

# FreeZone® Freeze Dryers & Accessories





## Table of Contents

|                                             |    |                                     |    |
|---------------------------------------------|----|-------------------------------------|----|
| <b>Why Labconco?</b>                        | 3  | <b>Accessories</b>                  |    |
| <b>Freeze Dryers At-a-Glance</b>            |    | Vacuum Pumps Oil & Filters          | 22 |
| Steps to Selecting a Freeze Dryer           | 4  | Product Shelves & Test Tube Racks   | 22 |
| FreeZone® Product Lines                     | 5  | Chambers                            | 23 |
| Lyo-Works™ Operating System                 |    | Heated Shelf Chambers               | 24 |
| Features and Benefits                       | 6  | Tray Dryer and Shelf Accessories    | 25 |
| <b>FreeZone® Freeze Dryers</b>              |    | Manifolds                           | 26 |
| Benchtop Freeze Dryers Features & Benefits  | 8  | End-Zone™ Starter Kit & Add-On Kit  | 27 |
| Benchtop Freeze Dryers Ordering Information | 9  | Stoppering Chambers & Optional Lids | 28 |
| Console Freeze Dryers Features & Benefits   | 10 | Secondary Traps                     | 29 |
| Console Freeze Dryers Ordering Information  | 11 | Fast-Freeze™ Flasks                 | 30 |
| Triad™ Freeze Dryers Features & Benefits    | 12 | Fast-Freeze™ Adapters               | 31 |
| Triad™ Freeze Dryers Ordering Information   | 13 | Fast-Freeze™ Accessories            | 32 |
| <b>FreeZone® Stoppering Tray Dryers</b>     |    | Serum Bottles & Ampules             | 33 |
| Features & Benefits                         | 14 | FreeZone® Benchtop Shell Freezers   | 34 |
| Standalone Stoppering Tray Dryers           | 15 | Carts, Shelves & Gauges             | 35 |
| Consoles w/Stoppering Tray Dryers           | 16 | <b>Specifications</b>               |    |
| <b>FreeZone® Bulk Tray Dryers</b>           | 17 | Freeze Dryers                       | 36 |
| <b>FreeZone® Small Tray Dryers</b>          | 18 | Triad™ Freeze Dryers                | 37 |
| <b>FreeZone® Lyph-Seal Tray Dryers</b>      | 19 | Stoppering Tray Dryers              | 38 |
| <b>Vacuum Pumps</b>                         | 20 | Bulk & Small Tray Dryers            | 39 |

# Built to last.

## The best materials.

Our mission is to build your equipment to the highest possible standards with superior quality materials.

- Environmentally responsible HCFC and CFC-free refrigerants with low GWP
- Durable powder-coated and stainless steel construction
- PTFE-coated collector chambers to resist corrosion
- Wear-resistant, quick-disconnect drain hoses

## The strongest pump protection.

Pump damage can be expensive and slow you down. That's why we make pump protection such a priority.

- Auto shut down if 5 mBar isn't reached during start up
- Patented\* drain line sensor detects moisture before it can damage your pump
- Vacuum break valve prevents oil backstreaming during power outages


## The most flexibility.

- Extend the life of your freeze dryer with endless accessories to help you adapt to the changing needs of your lab. No other manufacturer offers more
- Lyo-Works OS offers truly user-friendly touch screen controls, sample data monitoring and program storage. Stay updated every step of the way, in or out of the lab, with email alerts and our optional End-Zone end point detection



# FreeZone® Freeze Dryers

## Steps to Selecting a Freeze Dryer

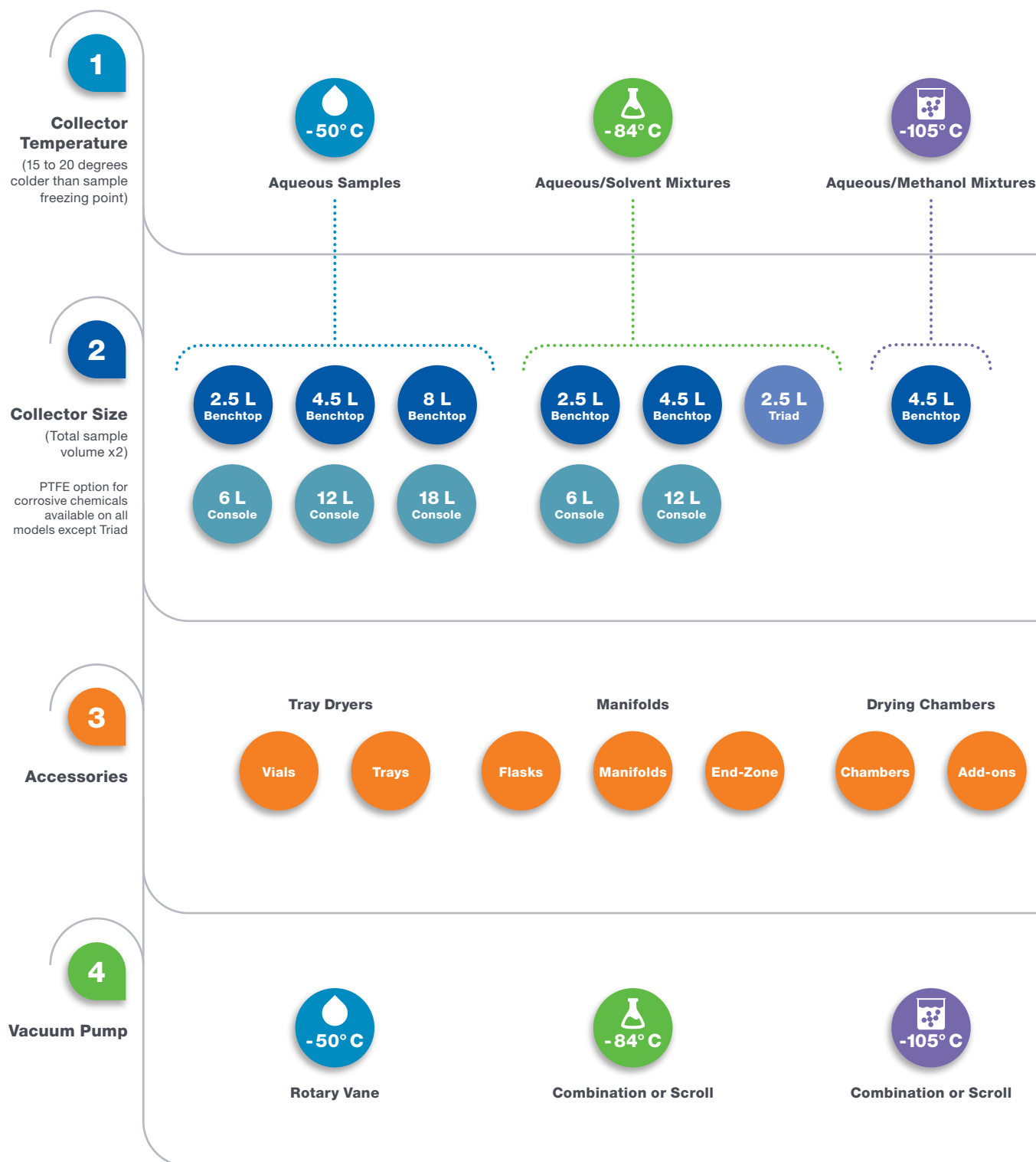

# FreeZone® Freeze Dryers

## FreeZone Product Lines

### Freeze Dryers



2.5 L

4.5 L

8 L

Benchtop



2.5 L

Triad Benchtop



6 L

12 L

18 L

Console and  
Console with Stoppering Tray Dryer

### Accessories




Tray Dryers



Manifolds



Drying Chambers

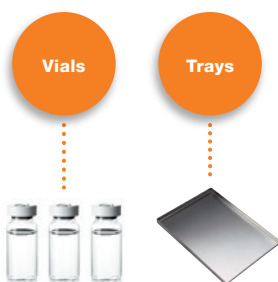

Vials

Trays



Flasks

Manifolds

End-Zone



Chambers

Add-ons

### Vacuum Pumps



Rotary Vane Pump



Chemically Resistant  
Rotary Vane Pump



Combination Pump



Scroll Pump



# Labconco for life.

Your research protects and improves how long, and how well, we live. To get there, you need confidence in the viability and stability of your samples.

At Labconco, we leave nothing to chance. You can expect our freeze dryers to deliver the absolute best quality, durability and flexibility available today.

You can also expect us to partner with you for the life of your freeze drying equipment, not just the warranty period. Which can be a long time, because we don't take shortcuts in design or production, and we don't use cut-rate materials.

Everyone talks about having better products and service, but with Labconco, it's actually true.

Labconco, for life.



# Lyo-Works™ Operating System

## Features & Benefits

### Data Logging

Starts/stops logging. View the data in table or graph form.

### Settings

Select language, security code, sensor calibration, time/date/time zone, alert emails and other settings.

### Sensors

Displays temperature values for sensors.

### Tray Dryers

Allows you to program or manually control shelf temperatures of installed tray dryers with temperature control.

### Programs

Lets you create, edit, view and initiate programs. Stores up to 30 programs with up to 16 steps each.

### Home Screen

Get at-a-glance collector temperature, vacuum pressure, and, if equipped, built-in shell freezer temperature.

### Auto

Initiates the start up sequence: collector cool down, vacuum pull down and optional purge valve. Ensures collector reaches -40° C before vacuum is applied.

### Collector

Manually starts/stops the collector and starts/stops defrost.

### Vacuum

Manually starts/stops the vacuum pump.

### End Point Detection

When End-Zone is installed, indicates when end point is reached for each sensor.

### Help


Access diagnostic tools, user manuals, videos, and updates

### Shell Freezer

Controls optional built-in shell freezer.

### Purge

Controls the optional purge valve.

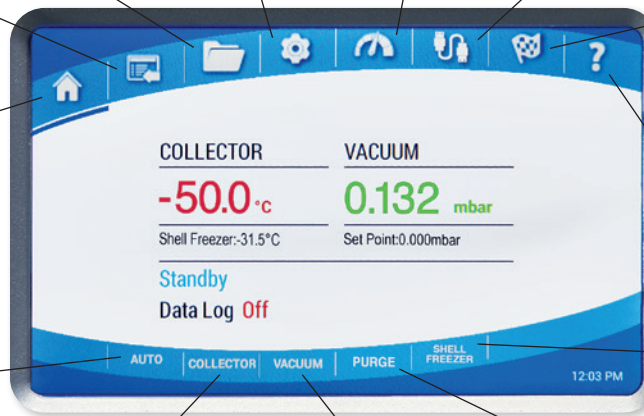

**Lyo-Works** is more than a large touch screen display. It's a set of powerful controls and useful insights about your freeze drying process. Lyo-Works gives you greater peace of mind about the integrity of your samples throughout lyophilization, even when you're not in the lab.



### USB port and Ethernet connection

Lets you easily output Lyo-Works data in CSV format to a flash drive or via email.

See how Lyo-Works works at [labconco.com/lyo-works](http://labconco.com/lyo-works).

# FreeZone® Benchtop Freeze Dryers

## Features & Benefits



### **Vacuum break valve**

prevents oil backstreaming when power to the freeze dryer or pump is shut off.

### **Patented\* moisture sensor**

protects the vacuum pump by preventing refrigeration or vacuum start-up when moisture is detected in collector chamber area.

### **Durable steel cabinet**

Benchtop models fit countertops and lab carts.

**USB port and Ethernet** to download onboard data.

### **Upright, stainless steel collector chamber**

speeds and simplifies defrost. Available with PTFE-coated collector coil and chamber for additional corrosion resistance.

### **Quick disconnect**

**drain hose fitting** for convenient defrost.

### **Vacuum control valve**

maintains set point vacuum level to speed the freeze drying process.

### **Hot gas defrost**

automatically shuts off when defrosting is complete.

**Full one year warranty**

\*U.S. Patent No. 6,968,631



# FreeZone® Benchtop Freeze Dryers

## Ordering Information



**2.5 Liter**

Light sample loads



**4.5 Liter**

Light to moderate sample loads



**8 Liter**

Moderate sample loads

**Catalog Number Configurator** Use this key to configure the nine digit catalog number.

7    0  0

Select the **REFRIGERATION TEMPERATURE** of your freeze dryer's collector. This is the second digit.

**0** = -50°C (-58°F). For samples containing water.

**1** = -84°C (-119°F). For low eutectic point samples containing water or acetonitrile.

**2** = -105°C (-157°F) at 60 Hz or -102°C (-152°F) at 50 Hz. For low eutectic point samples containing water, acetonitrile or dilute methanol or ethanol. (4.5 L only)

Select the **CAPACITY** of your freeze dryer's collector. These are the third and fourth digits.

**02** = 2.5 L      **04** = 4.5 L      **08** = 8 L (-50°C only)

Select the **COLLECTOR COIL**, **CHAMBER MATERIAL**, and **REFRIGERANT** of your freeze dryer. This is the sixth digit.

**1** = 1 = Stainless steel, Standard refrigerant

(Not available on 2.5L -50°C or 4.5L -50°C units)

**2** = PTFE-coated stainless steel, Standard refrigerant

(Not available on 2.5L -50°C or 4.5L -50°C units)

**5** = Stainless steel, Natural refrigerant

(Only available on 2.5L -50°C, 4.5L -50°C, and units ending in -90)

**6** = PTFE-coated stainless steel, Natural refrigerant

(Only available on 2.5L -50°C, 4.5L -50°C, and units ending in -90)

Select the **ELECTRICAL CONFIGURATION** and **OUTLET TYPE** of your freeze dryer. These are the eighth and ninth digits.

**00** = 115V, 60 Hz, with 20A plug for North America

**10** = 230V, 60 Hz, with plug for North America

**15** = 230V, 60 Hz, with plug for Saudi Arabia

**30** = 230V, 50 Hz, with Schuko plug

**40** = 230V, 50 Hz, with plug for China/Australia

**50** = 230V, 50 Hz, with plug for Great Britain (UK)

**70** = 230V, 50 Hz, with plug for India

**90** = 230V, 50 Hz, with Schuko, China/Australia, Great Britain (UK) and India

All models require a vacuum pump with proper displacement and 3/4" hose fitting, drying accessory, and freeze dry glassware if not bulk freeze drying

## Electrical Requirements, Dimensions & Shipping Weights

| Collector Capacity | Water removed in 24 hrs | Refrigeration Temperature | Volts | Amps | Overall Dimensions (w x d x h) | Shipping Weight  |
|--------------------|-------------------------|---------------------------|-------|------|--------------------------------|------------------|
| <b>2.5 L*</b>      | 2 L                     | -50°C (-58°F)             | 115   | 14** | 14.5" x 19.2" x 16.9"          | 120 lbs. (54 kg) |
| <b>2.5 L</b>       | 2 L                     | -50°C (-58°F)             | 230   | 7†   | 14.5" x 19.2" x 16.9"          | 120 lbs. (54 kg) |
| <b>2.5 L*</b>      | 2.2 L                   | -84°C (-119°F)            | 115   | 16** | 14.5" x 24.5" x 16.9"          | 165 lbs. (75 kg) |
| <b>2.5 L</b>       | 2.2 L                   | -84°C (-119°F)            | 230   | 10†  | 14.5" x 24.5" x 16.9"          | 165 lbs. (75 kg) |
| <b>4.5 L*</b>      | 2 L                     | -50°C (-58°F)             | 115   | 14** | 16.5" x 22.4" x 18.9"          | 139 lbs. (63 kg) |
| <b>4.5 L</b>       | 2 L                     | -50°C (-58°F)             | 230   | 7†   | 16.5" x 22.4" x 18.9"          | 139 lbs. (63 kg) |
| <b>4.5 L*</b>      | 4 L                     | -84°C (-119°F)            | 115   | 16** | 22.5" x 26.5" x 18.5"          | 197 lbs. (89 kg) |
| <b>4.5 L</b>       | 4 L                     | -84°C (-119°F)            | 230   | 10†  | 22.5" x 26.5" x 18.5"          | 197 lbs. (89 kg) |
| <b>4.5 L*</b>      | 2.5 L                   | -105°C (-15°F)            | 115   | 16** | 22.5" x 26.5" x 18.5"          | 205 lbs. (93 kg) |
| <b>4.5 L</b>       | 2.5 L                   | -105°C (-15°F)            | 230   | 10†  | 22.5" x 26.5" x 18.5"          | 205 lbs. (93 kg) |
| <b>8 L*</b>        | 4 L                     | -50°C (-58°F)             | 115   | 14*  | 22.5" x 26.5" x 18.5"          | 173 lbs. (79 kg) |
| <b>8 L</b>         | 4 L                     | -50°C (-58°F)             | 230   | 7**  | 22.5" x 26.5" x 18.5"          | 173 lbs. (79 kg) |

\*Requires 20 amp dedicated circuit.

\*\*System amperage shown includes 8 amp maximum vacuum pump rating. †System amperage shown includes 4.5 amp maximum vacuum pump rating.

# FreeZone® Console Freeze Dryers

## Features & Benefits



### Optional built-in shell freezer

with stainless steel bath and separate refrigeration system holds flasks up to 1200 mL. 6L only.

### Optional built-in mini vacuum drying chamber

holds small samples, either in bulk or in small containers such as serum bottles.

### Vacuum break valve

prevents oil backstreaming when power to the freeze dryer or pump is shut off.

### Durable steel cabinet

has casters for mobility and a two-piece brushed stainless steel front panel. Left side panel is easy to remove for access to the vacuum pump (pump sold separately). A window lets you see pump oil status.

### Patented\* moisture sensor

protects the vacuum pump by preventing refrigeration or vacuum start-up when moisture is detected in collector chamber area.

### Upright, stainless steel collector chamber

speeds and simplifies defrost. Available with PTFE-coated collector coil and chamber for additional corrosion resistance.

### Quick-disconnect drain hose fitting

for convenient defrost.

### Vacuum control valve

maintains set point vacuum level to speed the freeze drying process.

### Hot gas defrost

automatically shuts off when defrosting is complete.

### USB port and Ethernet

to download onboard data.

### Optional purge valve

isolates vacuum pump from the freeze dryer. The pump starts and warms the oil while the collector cools, before pulling vacuum. Can be used after run to purge contaminants from oil to extend pump life.

### Full one year warranty

\*U.S. Patent No. 6,968,631

# FreeZone® Console Freeze Dryers

## Ordering Information



6 Liter



12 Liter



18 Liter

- Moderate to large sample loads or numerous sample batches
- Compatible with the FreeZone Stoppering Tray Dryer or Bulk Tray Dryer (tray dryer sold separately)

**Catalog Number Configurator** Use this key to configure the nine digit catalog number.

7    1

Select the **REFRIGERATION TEMPERATURE** of your freeze dryer's collector. This is the second digit.

**0** = -50° C (-58° F). For samples containing water.

**1** = -84° C (-119° F). For low eutectic point samples containing water or acetonitrile. (6 L and 12 L only)

Select the **CAPACITY** of your freeze dryer's collector. These are the third and fourth digits.

**06** = 6 L      **12** = 12 L      **18** = 18 L (-50° C only)

Select the **COLLECTOR COIL**, **CHAMBER MATERIAL** and **REFRIGERANT** of your freeze dryer. Natural refrigerant only available on multi-plug electrical configurations (ending in -90). This is the sixth digit.

**1** = Stainless steel, Standard refrigerant

**2** = PTFE-coated stainless steel, Standard refrigerant

**5** = Stainless steel, Natural refrigerant

**6** = PTFE-coated stainless steel, Natural refrigerant

Select the **BUILT-IN OPTION**. This is the seventh digit.

**0** = None      **2** = Purge Valve and Mini Chamber.

**1** = Purge Valve      Add 20 lbs. (9 kg) to shipping weight.

**3** = Purge Valve and Shell Freezer. (6 L - 50° C only). Add 80 lbs. (36 kg) to shipping weight.

Select the **ELECTRICAL CONFIGURATION** and **OUTLET TYPE** of your freeze dryer. These are the eighth and ninth digits.

**00** = 115V, 60 Hz, with 20A plug for North America (6 L only)

**10** = 230V, 60 Hz, with plug for North America

**15** = 230V, 60 Hz, with plug for Saudi Arabia

**30** = 230V, 50 Hz, with Schuko plug

**40** = 230V, 50 Hz, with plug for China/Australia

**50** = 230V, 50 Hz, with plug for Great Britain (UK)

**70** = 230V, 50 Hz, with plug for India

**90** = 230V, 50 Hz, with Schuko, China/Australia, Great Britain (UK) and India

All models require a vacuum pump with proper displacement and 3/4" hose fitting, drying accessory, and freeze dry glassware if not bulk freeze drying.

## Electrical Requirements, Dimensions & Shipping Weights

| Collector Capacity | Water removed in 24 hrs | Refrigeration Temperature | Volts | Amps | Overall Dimensions (w x d x h) | Shipping Weight   |
|--------------------|-------------------------|---------------------------|-------|------|--------------------------------|-------------------|
| <b>6 L*</b>        | 4 L                     | -50° C (-58° F)           | 115   | 16** | 31.8" x 29.6" x 36.6"          | 285 lbs. (129 kg) |
| <b>6 L</b>         | 4 L                     | -50° C (-58° F)           | 230   | 8†   | 31.8" x 29.6" x 36.6"          | 285 lbs. (129 kg) |
| <b>6 L*</b>        | 4 L                     | -84° C (-119° F)          | 115   | 16** | 31.8" x 29.6" x 36.6"          | 299 lbs. (136 kg) |
| <b>6 L</b>         | 4 L                     | -84° C (-119° F)          | 230   | 10** | 31.8" x 29.6" x 36.6"          | 299 lbs. (136 kg) |
| <b>12 L</b>        | 8 L                     | -50° C (-58° F)           | 230   | 9†   | 31.8" x 29.6" x 36.6"          | 310 lbs. (141 kg) |
| <b>12 L</b>        | 6.5 L                   | -84° C (-119° F)          | 230   | 9†   | 31.8" x 29.6" x 36.6"          | 350 lbs. (159 kg) |
| <b>18 L</b>        | 10 L                    | -50° C (-58° F)           | 230   | 9*   | 31.8" x 29.6" x 36.6"          | 315 lbs. (143 kg) |

\*Requires 20 amp dedicated circuit.

\*\*System amperage shown includes 8 amp maximum vacuum pump rating. †System amperage shown includes 4.5 amp maximum vacuum pump rating.

# FreeZone® Triad® Freeze Dryers

## Features & Benefits



### Heated and cooled shelf

Separate HCFC/CFC-free refrigeration system and 1000 watt heater

### Four valve ports

Versatile Triad lets you freeze dry samples in flasks and on the shelf at the same time.

**USB port and Ethernet** to download onboard data.

### Pneumatic mechanism stoppers containers on shelf

Containers are stoppered while chamber is under vacuum and without the use of compressed gas. Atmospheric pressure causes the diaphragm to expand.

### Pre-freeze samples

on the shelves at a controlled rate (annealing). Improves sample quality and eliminates need for a separate freezer and product transfer.

### Large processing shelf

12.4" w x 14.5" d.

Accommodates containers with stoppers up to 148 mm maximum height.

### Collector drain pan and hose

Removeable stainless steel drain pan catches defrosted condensate.

**Full one year warranty**

## Serum Bottle Capacity

| Size   | Catalog Number | Shelf Capacity |
|--------|----------------|----------------|
| 2 ml   | <b>7575010</b> | 400            |
| 3 ml   | <b>7575210</b> | 441            |
| 5 ml   | <b>7573010</b> | 233            |
| 10 ml  | <b>7573210</b> | 196            |
| 20 ml  | <b>7573410</b> | 121            |
| 30 ml  | <b>7573610</b> | 86             |
| 50 ml  | <b>7573810</b> | 64             |
| 100 ml | <b>7574010</b> | 42             |
| 125 ml | <b>7574210</b> | 36             |

# FreeZone® Triad® Freeze Dryers

## Ordering Information



2.5 Liter Triad

- Small sample loads or numerous sample batches containing water or acetonitrile
- Stoppering vials or serum bottles under vacuum
- Limited floor space

**Catalog Number Configurator** Use this key to configure the nine digit catalog number.

7 9 4 0 0   0    

Select the **REFRIGERANT** of your freeze dryer. Natural refrigerant only available on multi-plug electrical configurations (ending in -90). This is the sixth digit.

**1** = Standard refrigerant      **5** = Natural refrigerant

Select the **ELECTRICAL CONFIGURATION** and **OUTLET TYPE** of your freeze dryer. These are the eighth and ninth digits.

**10** = 230V, 60 Hz, with plug for North America

**15** = 230V, 60 Hz, with plug for Saudi Arabia

**30** = 230V, 50 Hz, with Schuko plug

**40** = 230V, 50 Hz, with plug for China/Australia

**50** = 230V, 50 Hz, with plug for Great Britain (UK)

**70** = 230V, 50 Hz, with plug for India

**90** = 230V, 50 Hz, with Schuko, China/Australia, Great Britain (UK) and India

### All models require (not included):

- Vacuum pump with a displacement of at least 98 L/min, 0.02 mBar ultimate pressure and fitting for 3/4" ID vacuum hose. See pages 19-20
- Freeze dry glassware if not bulk freeze drying. See pages 29-32.
- Optional trays available. see page 21

### Dimensions & Shipping Weights

| Collector Capacity | Refrigeration Temperature               | Overall Dimensions (w x d x h) | Shipping Weight   |
|--------------------|-----------------------------------------|--------------------------------|-------------------|
| <b>2.5 L</b>       | -84° C (-119° F) Volts (230), Amps (12) | 28.3" x 29.5" x 27.7"          | 440 lbs. (200 kg) |

# FreeZone® Stoppering Tray Dryers

## Features & Benefits

### Heated and cooled shelves

Separate HCFC/CFC-free refrigeration system and 1000 watt heater.

### Vacuum release valve control

vents the chamber allowing the door to be opened.

**Stoppering knob** controls shelf stoppering mechanism.

**Three RTD sensor probes** monitor sample temperature.

### 1/8" OD backfill port

introduces inert gas from an outside source to the chamber, protecting samples from atmospheric moisture and contaminants.

### Pre-freeze samples

on the shelves at a controlled rate (annealing). Improves sample quality and eliminates need for a separate freezer and product transfer.

### Three large adjustable shelves

accommodate bulk trays or batches of serum bottles or ampules. Samples may be lyophilized and then stoppered on the shelves.



Shown with optional built-in 6-Port Tray Dryer Manifold

### Optional built-in 6-Port Tray Dryer Manifold




increases flexibility and capacity. Samples may be simultaneously freeze dried in flasks or other glassware. Also available as a field-installable accessory. See page 21 for ordering information.

### Optional built-in isolation valve

separates Stoppering Tray Dryer from Freeze Dryer to shorten time necessary to reestablish working vacuum levels and to provide a means for checking end point. Also available as a field-installable accessory. See page 21 for ordering information.

**Compatible with 6, 12, and 18L Console with Lyoworks OS.**

| Size   | Catalog Number | Shelf Capacity | Number of Shelves | Total Capacity |
|--------|----------------|----------------|-------------------|----------------|
| 2 ml   | <b>7575010</b> | 430            | 3                 | 1290           |
| 3 ml   | <b>7575210</b> | 483            | 3                 | 1449           |
| 5 ml   | <b>7573010</b> | 228            | 3                 | 684            |
| 10 ml  | <b>7573210</b> | 195            | 3                 | 585            |
| 20 ml  | <b>7573410</b> | 132            | 3                 | 396            |
| 30 ml  | <b>7573610</b> | 86             | 3                 | 258            |
| 50 ml  | <b>7573810</b> | 72             | 2                 | 144            |
| 100 ml | <b>7574010</b> | 42             | 2                 | 84             |
| 125 ml | <b>7574210</b> | 36             | 1                 | 36             |

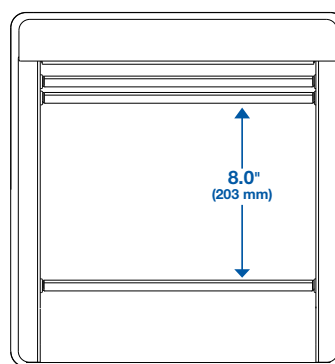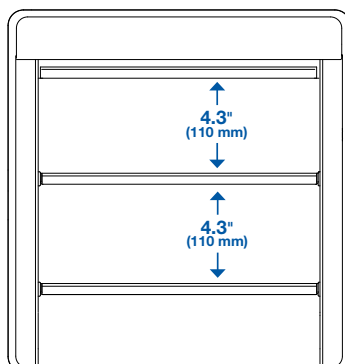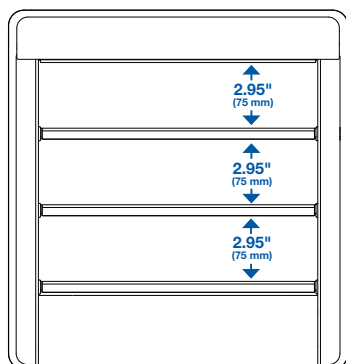

Stoppering Tray Dryer Shelf Spacing Dimensions.

Note: Adjustable shelves allow for additional spacing options.



# FreeZone® Stoppering Tray Dryers

## Ordering Information



Shown with optional built-in 6-Port Tray Dryer Manifold

**Catalog Number Configurator** Use this key to configure the nine digit catalog number.

To order your FreeZone Stoppering Tray Dryer and Freeze Dryer together as one catalog number, see page 16.

7 9 4 8 0        

Select **BUILT-IN OPTIONS**. This is the seventh digit.

**0** = None

**1** = 6-Port Tray Dryer Manifold

**2** = Isolation Valve

**3** = 6-Port Tray Dryer Manifold and Isolation Valve

Select the **REFRIGERANT** of your Stoppering Tray Dryer. Natural refrigerant only available on multi-plug electrical configurations (ending in -90). This is the sixth digit.

**1** = Standard refrigerant

**5** = Natural refrigerant

Select the **ELECTRICAL CONFIGURATION** and **OUTLET TYPE** of your freeze dryer. These are the eighth and ninth digits.

**00** =  115V, 60 Hz, with 20A plug for North America

**10** =  230V, 60 Hz, with plug for North America

**15** =  230V, 60 Hz, with plug for Saudi Arabia

**30** =  230V, 50 Hz, with Schuko plug

**40** =  230V, 50 Hz, with plug for China/Australia

**50** =  230V, 50 Hz, with plug for Great Britain (UK)

**70** =  230V, 50 Hz, with plug for India

**90** =     230V, 50 Hz, with Schuko, China/Australia, Great Britain (UK) and India

### All models require (not included):

- FreeZone 6, 12 or 18 Liter Console Freeze Dryer with Lyo-Works OS. FreeZone Stoppering Tray Dryer and Freeze Dryer may be ordered as separate components or together as one catalog number on selected models of FreeZone 6, 12 or 18 Liter Console Freeze Dryers. See page 16
- Vacuum pump with a displacement of at least 163 L/min, 0.02 mBar ultimate pressure and fitting for 3/4" ID vacuum hose. See pages 19-20
- Freeze dry glassware if not bulk freeze drying. See pages 29-32

### Electrical Requirements

| Volts | Amps |
|-------|------|
| 115   | 16   |
| 230   | 9    |

### Shipping Weights

| Built-In Options Included                      | Shipping Weight   |
|------------------------------------------------|-------------------|
| None                                           | 423 lbs. (192 kg) |
| 6-Port Tray Dryer Manifold                     | 438 lbs. (199 kg) |
| Isolation Valve                                | 431 lbs. (196 kg) |
| 6-Port Tray Dryer Manifold and Isolation Valve | 458 lbs. (208 kg) |

# FreeZone® Console Freeze Dryers with Stoppering Tray Dryer

## Ordering Information



6 Liter



12 Liter



18 Liter

- Moderate to large sample loads or numerous sample batches
- Stoppering vials or serum bottles under vacuum

**Catalog Number Configurator** Use this key to configure the nine digit catalog number.

7    2

Select the **REFRIGERATION TEMPERATURE** of your freeze dryer's collector. This is the second digit.

**0** = -50° C (-58° F). For samples containing water.

**1** = -84° C (-119° F). For low eutectic point samples containing water or acetonitrile. (6 L and 12 L only)

Select the **CAPACITY** of your freeze dryer's collector. These are the third and fourth digits.

**06** = 6 L      **12** = 12 L      **18** = 18 L (-50° C only)

Select the **COLLECTOR COIL, CHAMBER MATERIAL, and REFRIGERANT** of your freeze dryer. Natural refrigerant only available on multi-plug electrical configurations (ending in -90). This is the sixth digit.

**1** = Stainless steel, Standard refrigerant

**2** = PTFE-coated stainless steel, Standard refrigerant

**5** = Stainless steel, Natural refrigerant

**6** = PTFE-coated stainless steel, Natural refrigerant

Select the **BUILT-IN OPTION**. This is the seventh digit.

**0** = None      **1** = Purge Valve

Select the **ELECTRICAL CONFIGURATION** and **OUTLET TYPE** of your freeze dryer. These are the eighth and ninth digits.

**00** = 115V, 60 Hz, with 20A plug for North America (6 L)

**10** = 230V, 60 Hz, with plug for North America

**15** = 230V, 60 Hz, with plug for Saudi Arabia

**30** = 230V, 50 Hz, with Schuko plug

**40** = 230V, 50 Hz, with plug for China/Australia

**50** = 230V, 50 Hz, with plug for Great Britain (UK)

**70** = 230V, 50 Hz, with plug for India

**90** = 230V, 50 Hz, with Schuko, China/Australia, Great Britain (UK) and India

All models require a vacuum pump with proper displacement and 3/4" hose fitting and freeze dry glassware if not bulk freeze drying.

## Electrical Requirements, Dimensions (Tray Dryer mounted on Console) & Shipping Weights (Two Cartons)

| Collector Capacity | Water removed in 24 hrs | Refrigeration Temperature | Volts | Console Amps | Tray Dryer Amps | Overall Dimensions (w x d x h) | Shipping Weight   |
|--------------------|-------------------------|---------------------------|-------|--------------|-----------------|--------------------------------|-------------------|
| <b>6 L*</b>        | 4 L                     | -50° C                    | 115   | 16**         | 16              | 32.4" x 29.6" x 63.7"          | 678 lbs. (308 kg) |
| <b>6 L</b>         | 4 L                     | -50° C                    | 230   | 8†           | 9               | 32.4" x 29.6" x 63.7"          | 678 lbs. (308 kg) |
| <b>6 L*</b>        | 4 L                     | -84° C                    | 115   | 16**         | 16              | 32.4" x 29.6" x 63.7"          | 693 lbs. (314 kg) |
| <b>6 L</b>         | 4 L                     | -84° C                    | 230   | 10**         | 9               | 32.4" x 29.6" x 63.7"          | 693 lbs. (314 kg) |
| <b>12 L</b>        | 8 L                     | -50° C                    | 230   | 9†           | 9               | 32.4" x 29.6" x 63.7"          | 698 lbs. (317 kg) |
| <b>12 L</b>        | 6.5 L                   | -84° C                    | 230   | 9†           | 9               | 32.4" x 29.6" x 63.7"          | 743 lbs. (337 kg) |
| <b>18 L</b>        | 10 L                    | -50° C                    | 230   | 9†           | 9               | 32.4" x 29.6" x 63.7"          | 713 lbs. (323 kg) |

\*Requires 20 amp dedicated circuit.

\*\*System amperage shown includes 8 amp maximum vacuum pump rating. †System amperage shown includes 4.5 amp maximum vacuum pump rating.

# FreeZone® Bulk Tray Dryers

## Ordering Information



### 3 Heated Shelves with Sensors

Each stainless steel shelf includes a 200-watt heater and a sample temperature probe.

### Bulk Tray Dryers with heated shelves

Moderate to large bulk sample loads or numerous sample batches

### Bulk Tray Dryer models require (not included):

- FreeZone 6, 12 or 18 Liter Console Freeze Dryer with Lyo-Works OS. See pages 10-11
- Vacuum pump with a displacement of at least 163 L/min, 0.02 mBar ultimate pressure and fitting for 3/4" ID vacuum hose. See pages 19-20
- Freeze dry glassware if not bulk freeze drying. See pages 29-32
- 115V models require 8A; 230V models require 4A

**Catalog Number Configurator** Use this key to configure the nine digit catalog number.

7 8 0 6 0 1      

Select **BUILT-IN OPTIONS**. This is the seventh digit.

- 0 = None      1 = 6-Port Tray Dryer Manifold  
2 = Isolation Valve      3 = 6-Port Tray Dryer Manifold & Isolation Valve

Select the **ELECTRICAL CONFIGURATION** and **OUTLET TYPE**. These are the eighth and ninth digits.

- 00 = 115V, 60 Hz, with 20A plug for North America  
10 = 230V, 60 Hz, with plug for North America  
15 = 230V, 60 Hz, with plug for Saudi Arabia  
30 = 230V, 50 Hz, with Schuko plug  
40 = 230V, 50 Hz, with plug for China/Australia  
50 = 230V, 50 Hz, with plug for Great Britain (UK)  
70 = 230V, 50 Hz, with plug for India  
90 = 230V, 50 Hz, with Schuko, China/Australia, Great Britain (UK) and India


## Shipping Weights

| Options Included                               | Shipping Weight   |
|------------------------------------------------|-------------------|
| None                                           | 248 lbs. (113 kg) |
| 6-Port Tray Dryer Manifold                     | 263 lbs. (119 kg) |
| Isolation Valve                                | 256 lbs. (116 kg) |
| 6-Port Tray Dryer Manifold and Isolation Valve | 283 lbs. (128 kg) |

### Optional Heated Shelves with Sensors 12.7" w x 16.6" d.

FreeZone Bulk Tray Dryers include three and can accommodate two additional Heated Shelves with Sensors.

| Catalog Number | Description                        | Shipping Weight   |
|----------------|------------------------------------|-------------------|
| 730151020      | Shelf for 115V, 50/60 Hz operation | 6.0 lbs. (2.7 kg) |
| 730151030      | Shelf for 230V, 50/60 Hz operation | 6.0 lbs. (2.7 kg) |

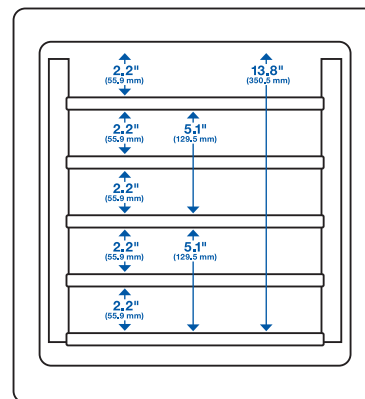

Bulk Tray Dryer Shelf Spacing Dimensions



# FreeZone® Small Tray Dryers

## Ordering Information



### Small Tray Dryers with heated shelves

Small to moderate bulk sample loads

### Heated Small Shelves with Sensors

Includes 3 Heated Shelves with option for 2 additional heated shelves. Each stainless steel shelf includes a 200-watt heater and a sample temperature probe


**Optional Heated Small Shelves with Sensors** 10.6" w x 12.0" d. FreeZone Small Tray Dryers include three and can accommodate two additional Heated Shelves with Sensors.

| Catalog Number   | Description                                                              | Shipping Weight  |
|------------------|--------------------------------------------------------------------------|------------------|
| <b>730152020</b> | Shelf for 115V, 50/60 Hz operation                                       | 15 lbs. (7 kg)   |
| <b>730152030</b> | Shelf for 230V, 50/60 Hz operation                                       | 15 lbs. (7 kg)   |
| <b>736581000</b> | Temperature Probe only<br>(For Bulk and Small Tray Dryer Heated Shelves) | 0.5 lb. (0.2 kg) |

| Catalog Number   | Electrical Requirement | Receptacle & Plug Type                                                 |
|------------------|------------------------|------------------------------------------------------------------------|
| <b>780701000</b> | 115V, 60 Hz, 5A        | 115V, 60 Hz for North America                                          |
| <b>780701010</b> | 230V, 60 Hz, 2.5A      | 230V, 60 Hz for North America                                          |
| <b>780701015</b> | 230V, 60 Hz, 2.5A      | 230V, 60 Hz for Saudi Arabia                                           |
| <b>780701030</b> | 230V, 50 Hz, 2.5A      | 230V, 50 Hz Schuko                                                     |
| <b>780701040</b> | 230V, 50 Hz, 2.5A      | 230V, 50 Hz for China/Australia                                        |
| <b>780701050</b> | 230V, 50 Hz, 2.5A      | 230V, 50 Hz for Great Britain (UK)                                     |
| <b>780701070</b> | 230V, 50 Hz, 2.5A      | 230V, 50 Hz for India                                                  |
| <b>780701090</b> | 230V, 50 Hz, 2.5A      | 230V, 50 Hz for Schuko, China/Australia, Great Britain (UK), and India |

### Small Tray Dryer models require (not included):

- Any FreeZone Freeze Dryer with Lyo-Works OS
- Vacuum pump with a displacement of at least 98 L/min, 0.020 mBar ultimate pressure and fitting for 3/4" ID vacuum hose. See pages 20-21
- Freeze dry glassware if not bulk freeze drying. See pages 29-32

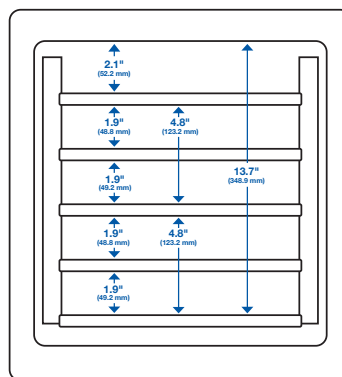

**7522700—Dry Ice Benchtop Freeze Dryer.** For occasional lyophilization requirements. Chamber, 8.8" dia. x 9.8" h type 304 stainless steel with 12 valves and capacity of 1 liter of ice. Valves accommodate either 1/2" or 3/4" flask adapters. Includes dry ice/solvent center well with 1.92 liter capacity and cover, 1/2" OD port for connection to vacuum pump and 3/8" OD port for connection to vacuum gauge. **Dry ice, heat transfer solution, vacuum pump, vacuum gauge, tubing and glassware are required (not included).** See pages 19-20 for vacuum pumps and accessories. See pages 29-32 for glassware. Shipping weight 11 lbs. (5 kg)



# FreeZone™ Lyph-Seal Tray Dryers

## Ordering Information



### Versatile Sample Preservation

Protect your samples from degradation by freeze drying and sealing under vacuum in the same device. Bulk samples that are too large for vials or ampules can now be preserved within flexible, zipper seal bags.

### Temperature Controlled Shelves

Four independently heated sample shelves drive the lyophilization process just as you specify with temperature setpoints between -20° C and 60° C in 1-degree increments. Shelves are removable to accommodate larger samples.

### Lyph-Seal Tray Dryer models require (not included):

- FreeZone Freeze Dryer with Lyo-Works OS
- Vacuum pump with a displacement of at least 98 L/min, 0.020 mBar ultimate pressure and fitting for 3/4" ID vacuum hose. See pages 20-21
- Lyph-Seal Sample Bags

| Catalog Number | Electrical Requirement | Receptacle & Plug Type                                                 |
|----------------|------------------------|------------------------------------------------------------------------|
| 780401000      | 115V, 60 Hz, 4A        | 115V, 60 Hz for North America                                          |
| 780401010      | 230V, 60 Hz, 2A        | 230V, 60 Hz for North America                                          |
| 780401015      | 230V, 60 Hz, 2A        | 230V, 60 Hz for Saudi Arabia                                           |
| 780401030      | 230V, 50 Hz, 2A        | 230V, 50 Hz Schuko                                                     |
| 780401040      | 230V, 50 Hz, 2A        | 230V, 50 Hz for China/Australia                                        |
| 780401050      | 230V, 50 Hz, 2A        | 230V, 50 Hz for Great Britain (UK)                                     |
| 780401070      | 230V, 50 Hz, 2A        | 230V, 50 Hz for India                                                  |
| 780401090      | 230V, 50 Hz, 2A        | 230V, 50 Hz for Schuko, China/Australia, Great Britain (UK), and India |




### Lyph-Seal Tray Dryer Accessories

**7592100—Lyph-Seal Sample Bags** 6" W x 10.25" H  
Sanitary zipper seal bags made from durable PET/PE. Designed for use with the Lyph-Seal Tray Dryer system and compatible with commercially available impulse sealers. Pack of 100.

### 7594700—Lyph-Seal Bag Freezer Form

6.7" W x 7.7" D x 6.7" H

Stainless steel form ensures uniform sample thickness during pre-freezing.

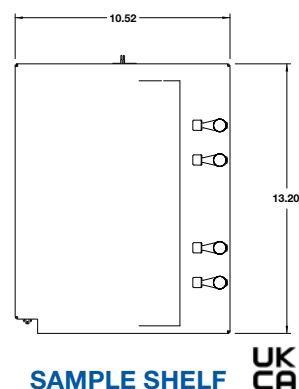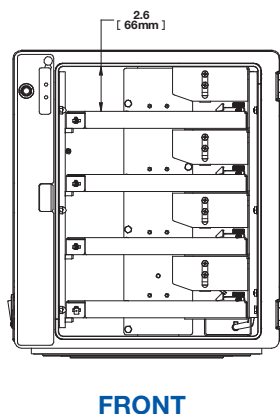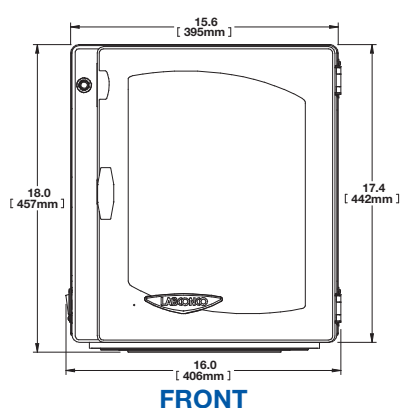


# Vacuum Pumps

## Features & Benefits



### Scroll Vacuum Pumps

- Designed for use with acids below 20% concentration and solvents including acetonitrile
- Two spiral scrolls compress air and vapors to move them toward the exhaust
- Oil-free design has no hydrocarbons
- Ultimate vacuum  $7 \times 10^{-3}$  mBar (5.3 micron) for 7587000
- Ultimate vacuum (partial pressure)  $2 \times 10^{-2}$  mBar (15 micron) for 7587100

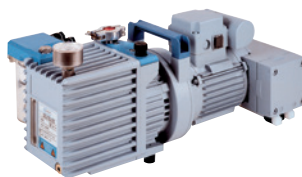

### VACUUBRAND HYBRID\* Vacuum Pumps

- Designed for use with acids and other harsh chemicals including TFA, TFA by-products, acetonitrile, HBe and  $\text{HNO}_3$  present in samples such as HPLC-prepared and peptide purified materials
- Comprised of a two-stage rotary vane pump and two-stage, chemical-resistant diaphragm pump. Superior solvent and acid handling capabilities
- Ultimate vacuum  $2 \times 10^{-3}$  mBar (1.5 micron)



### Chemically Resistant Rotary Vane Vacuum Pump

- Designed for use with harsh chemicals such as TFA, acetonitrile, HBr, and others
- System includes continuous acid neutralization and oil filtration, extending the service interval of the pump in caustic conditions
- Large oil chamber slows oil breakdown within pump
- Ultimate vacuum  $2 \times 10^{-3}$  mBar (1.5 micron)



### General Purpose Rotary Vane Vacuum Pumps

- Designed for use with aqueous samples
- Single phase direct drive motor, totally enclosed and fan cooled. Should the motor overheat, the thermal overload device switches off the pump. When the pump cools down, the motor automatically restarts
- Ultimate vacuum  $2 \times 10^{-3}$  mBar (1.5 micron)

| FreeZone Model                                                                      |                                                                                                     | Sample Contains  | Pump Type      | Recommended Pump<br>(≥98 lpm)** |           | Recommended Pump<br>(≥130 lpm) |           |
|-------------------------------------------------------------------------------------|-----------------------------------------------------------------------------------------------------|------------------|----------------|---------------------------------|-----------|--------------------------------|-----------|
|                                                                                     |                                                                                                     |                  |                | 115 Volts                       | 230 Volts | 115 Volts                      | 230 Volts |
|  | Benchtop 2.5, 4.5, & 8<br>Liter or<br>Dry Ice Freeze Dryer                                          | Water            | Rotary Vane    | 1472100                         | 7739402   | 1467700                        | 7739403   |
|                                                                                     |                                                                                                     | Acids & Solvents | CR Rotary Vane | 7584500                         | 7584501   | 7584500                        | 7584501   |
|                                                                                     |                                                                                                     | Acids & Solvents | Hybrid         | 7584000                         | 7584002   |                                |           |
|                                                                                     |                                                                                                     | Acids & Solvents | Scroll         | 7587100                         | 7587100   | 7587000                        | 7587000   |
|  | Console 6, 12 & 18 Liter                                                                            | Water            | Rotary Vane    | 1472100                         | 7739402   | 1467700                        | 7739403   |
|                                                                                     |                                                                                                     | Acids & Solvents | CR Rotary Vane | 7584500                         | 7584501   | 7584500                        | 7584501   |
|                                                                                     |                                                                                                     | Acids & Solvents | Hybrid         | 7584000                         | 7584002   |                                |           |
|                                                                                     |                                                                                                     | Acids & Solvents | Scroll         | 7587100                         | 7587100   | 7587000                        | 7587000   |
|  | Console 6, 12 & 18 Liter<br>with Stoppering Tray Dryer,<br>Bulk Tray Dryer<br>or Tall Clear Chamber | Water            | Rotary Vane    |                                 |           | 1467700                        | 7739403   |
|                                                                                     |                                                                                                     | Acids & Solvents | CR Rotary Vane |                                 |           | 7584500                        | 7584501   |
|                                                                                     |                                                                                                     | Acids & Solvents | Scroll         |                                 |           | 7587000                        | 7587000   |
|  | Triad                                                                                               | Water            | Rotary Vane    |                                 | 7739402   |                                | 7739403   |
|                                                                                     |                                                                                                     | Acids & Solvents | CR Rotary Vane |                                 | 7584501   |                                | 7584501   |
|                                                                                     |                                                                                                     | Acids & Solvents | Hybrid         |                                 | 7584002   |                                |           |
|                                                                                     |                                                                                                     | Acids & Solvents | Scroll         |                                 | 7587100   |                                | 7587000   |

\*VACUUBRAND™ and VACUUBRAND HYBRID™ are trademarks of VACUUBRAND GMBH + CO KG.

\*\*Faster freeze drying rates may be achieved by using pumps with larger displacement capacities.


# Vacuum Pumps

## Ordering Information



### Scroll Vacuum Pumps

| Catalog Number | Electrical Configuration                                              | Displacement (L/min) | Dimensions (w x d x h) | Shipping Weight |
|----------------|-----------------------------------------------------------------------|----------------------|------------------------|-----------------|
| <b>7587000</b> | 115V/ 230V, 50/60 Hz, 10A/6A<br>(Includes NEMA 5-15P and reverse IEC) | 190                  | 20.0" x 12.0" x 15.0"  | 59 lbs. (27 kg) |
| <b>7587100</b> | 115V/ 230V, 50/60 Hz, 10A/6A<br>(Includes NEMA 5-15P and reverse IEC) | 103                  | 20.0" x 12.0" x 15.0"  | 59 lbs. (27 kg) |

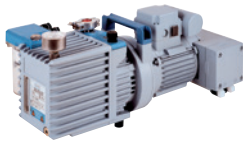

### VACUUBRAND HYBRID\*\* Vacuum Pumps

| Catalog Number | Electrical Configuration                         | Displacement at 60 Hz (50 Hz) L/min | Dimensions (w x d x h) | Shipping Weight |
|----------------|--------------------------------------------------|-------------------------------------|------------------------|-----------------|
| <b>7584000</b> | 115V, 60 Hz, 5.7A (Includes 115V NEMA 5-15P)     | 115 (98)                            | 20.0" x 12.0" x 9.1"   | 66 lbs. (30 kg) |
| <b>7584002</b> | 230V, 50/60 Hz, 2.7A (Includes 230V reverse IEC) | 115 (98)                            | 20.0" x 12.0" x 9.1"   | 66 lbs. (30 kg) |



### Chemically Resistant Rotary Vane Vacuum Pumps

| Catalog Number | Electrical Configuration                    | Displacement at 60 Hz (50 Hz) L/min | Dimensions (w x d x h) | Shipping Weight     |
|----------------|---------------------------------------------|-------------------------------------|------------------------|---------------------|
| <b>7584500</b> | 115V, 60Hz, 4.3A (includes 115V NEMA 5-15P) | 158                                 | 18.2" x 11.1" x 9.2"   | 61.6 lbs. (27.9 kg) |
| <b>7584501</b> | 230V, 50/60Hz, 2.1A (includes reverse IEC)  | 158 (133)                           | 18.2" x 11.1" x 9.2"   | 61.6 lbs. (27.9 kg) |



### General Purpose Rotary Vane Vacuum Pumps

| Catalog Number | Electrical Configuration                         | Displacement at 60 Hz (50 Hz) L/min | Dimensions (w x d x h) | Shipping Weight |
|----------------|--------------------------------------------------|-------------------------------------|------------------------|-----------------|
| <b>1472100</b> | 115V, 50/60 Hz, 4.6A (Includes 115V NEMA 5-15P)  | 117 (98)                            | 6.5" x 17.0" x 10.4"   | 62 lbs. (28 kg) |
| <b>7739402</b> | 230V, 50/60 Hz, 2.4A (Includes 230V reverse IEC) | 117 (98)                            | 6.5" x 17.0" x 10.4"   | 62 lbs. (28 kg) |
| <b>1467700</b> | 115V, 50/60 Hz, 7.8A (Includes 115V NEMA 5-15P)  | 195 (163)                           | 6.5" x 18.5" x 10.4"   | 65 lbs. (30 kg) |
| <b>7739403</b> | 230V, 50/60 Hz, 4.0A (Includes 230V reverse IEC) | 195 (163)                           | 6.5" x 18.5" x 10.4"   | 65 lbs. (30 kg) |

\*\*VACUUBRAND HYBRID™ is a trademark of VACUUBRAND GMBH + CO KG.

# Accessories

## Vacuum Pump Oil, Filters and Accessories

### Vacuum Pump Oil

| Cat. No. | Product Name         | Volume            | Shipping Weight  | Compatibility       |
|----------|----------------------|-------------------|------------------|---------------------|
| 1988000  | Vacuum Pump Oil      | One liter         | 2.0 lbs. (0.9kg) | Rotary Vane Pump    |
| 1988001  | Vacuum Pump Oil      | One gallon (3.8L) | 12 lbs. (5kg)    | Rotary Vane Pump    |
| 7772700  | Gold Vacuum Pump Oil | One liter         | 2.0 lbs. (0.9kg) | CR Rotary Vane Pump |
| 7584300  | VACUUBRAND B-Oil     | One liter         | 2.0 lbs. (0.9kg) | Hybrid Pump         |
| 7541300  | Vacuum Pump Oil      | One gallon (3.8L) | 12 lbs. (5kg)    | Hybrid Pump         |

**1473500—Oil Return Kit.** For Rotary Vane Pump. Allows extended operation of the pump with the ballast open and helps to keep oil purged of contaminants. Includes stainless steel cylinder and tubing with clamp. Shipping weight 0.5 lb. (0.2 kg)

### Vacuum Pump Filters

| Cat. No. | Product Name                                  | Shipping Weight  | Compatibility       |
|----------|-----------------------------------------------|------------------|---------------------|
| 1472200  | Inlet Filter Assembly and Cartridge           | 1.0 lb. (0.5kg)  | Rotary Vane Pump    |
| 1472500  | Replacement Inlet Filter Cartridge            | 0.5 lb. (0.2kg)  | Rotary Vane Pump    |
| 1473400  | Replacement Exhaust Filter                    | 1.0 lb. (0.5kg)  | Rotary Vane Pump    |
| 1473200  | Replacement Oil Mist Exhaust Filter Cartridge | 0.5 lb. (0.2 kg) | Rotary Vane Pump    |
| 1473300  | Replacement Odor Exhaust Filter Cartridge     | 0.5 lb. (0.2 kg) | Rotary Vane Pump    |
| 7584600  | Alumina Filter Element for Oil Sump           | 1.0 lb. (0.5kg)  | CR Rotary Vane Pump |

**7587300 - Pump Silencer.** Filters the exhaust while silencing exhaust pump noise. Includes filter, centering ring and clamp to attach to exhaust outlet. Shipping weight 2.0 lbs. (0.9 kg)

**7587200 - Tip Seal Kit for Scroll Vacuum Pump.** Tip Seal Replacement Kit provides all parts needed to replace tip seals in Scroll Vacuum Pumps. Includes two tip seals and O-rings. Labconco recommends replacing tip seals every 40,000 hours of pump use. Shipping weight 0.5 lbs (0.2 kg)

## Product Shelves & Test Tube Racks



**Product Shelves.** For bulk samples, microwell plates or serum bottles. Aluminum 10" dia. with 1/4" rims, three steel support rods and rubber feet. **Requires Chamber (not included).**

**7317100—Three.** Overall height: 9.8".  
Shipping weight 4.5 lbs. (2.0 kg)

**7317200—Five.** Overall height: 17.5".  
Shipping weight 7.0 lbs. (3.2 kg)



**Round Test & Microcentrifuge Tube Racks.** 304 stainless steel. 9.0" dia. x 2.5" h. Holds up to 104 each 12-13 mm test tubes or microcentrifuge tubes. For use with 7317100, 7317200 and Heated Product Shelf Chambers with nine-digit catalog numbers. **Tubes not included.**

**7379000—Round Test Tube Rack.** Shipping weight 2.8 lbs. (1.3 kg)

**7379001—Microcentrifuge Tube Rack.** Includes adapter plate for microcentrifuge tubes. Shipping weight 2.8 lbs. (1.3 kg)

# Accessories

## Chambers



**304 Stainless Steel Drying Chambers.** Clear lid with neoprene gasket and valves and plastic molded knobs. For 1/2" and 3/4" flask adapters.

**7522800—12-Port.** 304 stainless steel chamber with 1/2" thick acrylic lid and twelve neoprene valves. 8.0" dia. x 9.75" h. Shipping weight 11 lbs. (5 kg)

**7522801—12-Port.** PTFE-coated 304 stainless steel chamber with 1/2" thick acrylic lid and twelve neoprene valves. 8.0" dia. x 9.75" h. Shipping weight 11 lbs. (5 kg)

**7522900—16-Port.** 304 stainless steel chamber with 3/4" thick acrylic lid and sixteen neoprene valves. 13.0" dia. x 13.0" h. Shipping weight 21 lbs. (10 kg)

**7522901—16-Port.** PTFE-coated 304 stainless steel chamber with 3/4" thick acrylic lid and sixteen neoprene valves. 13.0" dia. x 13.0" h. Shipping weight 21 lbs. (10 kg)

**7555200—0-Port.** 304 stainless steel chamber with 3/4" thick glass lid, vacuum release valve, and pass-through for electrical cords. For use with 7317100 product shelves. 13.0" dia. x 13.0" h. Shipping weight 21 lbs. (10 kg)



**Clear Chambers with Valves.** For bulk drying with flask drying option. Includes acrylic base, anodized aluminum lid, and eight neoprene valves with molded plastic knobs. For 1/2" and 3/4" flask adapters. 8" dia. opening and cap for 3" dia. attachment port.

**7318802—Short.** For use with 7317100 product shelves. 12.0" ID dia. x 15.0" h. Shipping weight 22 lbs. (10 kg)

**7318800—Tall.** For use with 7317100 and 7317200 product shelves. 12.0" ID dia. x 22.8" h. Shipping weight 25 lbs. (11 kg)



**Clear Chambers.** For view of product shelves during bulk drying. Includes acrylic lid and base, neoprene gasket, vacuum release valves and pass-through for electrical cords. 8" dia. opening and cap for 3" dia. attachment port on console freeze dryer.

**7318702—Short.** For use with 7509200, 7509201 and 7317100 product shelves 12.0" ID dia. x 16.5" h. Shipping weight 17 lbs. (8 kg)

**7318700—Tall.** For use with 7509400, 7509401 and 7317200 product shelves 12.0" ID dia. x 24.3" h. Shipping weight 20 lbs. (10 kg)

# Drying Accessories

## Heated Shelf Chambers



**Heated Product Shelf Chambers with Valves.** For bulk drying and flask drying. Eight ports accommodate both 1/2" and 3/4" flask adapters. Includes 12" dia. acrylic cylinder and base, 10" dia. removable shelves, and temperature probes. Controlled heat settings -20 C (-4F) to + 60 C (+140F). Shelves do not cool. Overall dimensions: 13.5" w x 16.8" d x 14.2" h or 21.8" h. Has 8" dia. opening and cap for 3" dia. attachment port.



**Heated Product Shelf Chambers.** For bulk drying solid samples or drying small serum bottles or vials. Include 12" dia. acrylic cylinder and base, 10" dia. removable shelves, temperature probes and vacuum release valve. Controlled heat settings -20 C (-4F) to + 60 C (+140F). Shelves do not cool. Overall dimensions: 13.5" w x 15.2" d. x 14.2" h or 21.8" h. Has 8" dia. opening and cap for 3" dia. attachment port.

7 8 0 8 0        

Select of **NUMBER OF VALVES**. This is the sixth digit.

- 1** = Zero Valves
- 2** = 8 Valves

Select of **NUMBER OF SHELVES**. This is the seventh digit.


- 0** = 3 shelves
- 1** = 5 shelves

Select the **ELECTRICAL CONFIGURATION** and **OUTLET TYPE** of your freeze dryer. These are the eighth and ninth digits.

- 00** = 115V, 60 Hz, with 20A plug for North America
- 10** = 230V, 60 Hz, with plug for North America
- 15** = 230V, 60 Hz, with plug for Saudi Arabia
- 30** = 230V, 50 Hz, with Schuko plug
- 40** = 230V, 50 Hz, with plug for China/Australia
- 50** = 230V, 50 Hz, with plug for Great Britain (UK)
- 70** = 230V, 50 Hz, with plug for India
- 90** = 230V, 50 Hz, with Schuko, China/Australia, Great Britain (UK) and India

# Accessories

## Tray Dryer and Shelf

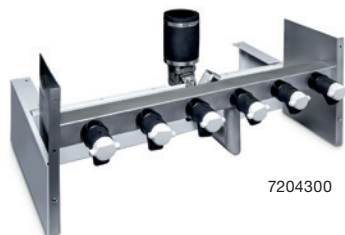


7204300

**7204300—6-Port Manifold with Isolation Valve.** For 1/2" or 3/4" flask adapters. **Not for use with Small Tray Dryers.** 30.1" w x 20.6" d x 13.5" h. Shipping weight 45 lbs. (21 kg)

**7726500—6-Port Tray Dryer Manifold.** Not field installable with 7761500. **Not for use with Small Tray Dryers.** 30.0" w x 16.5" d x 7.0" h. Shipping weight 19 lbs. (9 kg)

**7761500—Isolation Valve.** Separates tray and freeze dryers to shorten time necessary to reestablish working vacuum levels and provide a means for checking end point. Includes valve and handle, clamps and couplings. Not field installable with 7726500. **Not for use with Small Tray Dryers.** Shipping weight 5.0 lbs. (2.3 kg)

**7756300—Shelf Spacers.** Assist stoppering of partial loads of small serum bottles. 304 stainless steel. One pair. 9.0" w x 13.0" d x 2.5" h. Shipping weight 3.0 lbs. (1.4 kg)

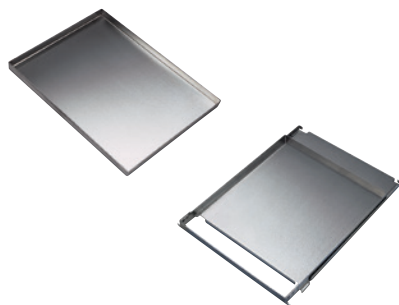

**Bulk Trays.** Contain pre-frozen liquids for bulk drying. 304 stainless steel with 0.9" lip. **Large Tray not for use with Small Tray Dryer.** Shipping weight 3.0 lbs. (1.4 kg)

**7377900—Small.** 10.3" x 11.6"

**7756200—Large.** 12.4" x 14.5" Compatible with Stoppering Tray Dryer, Bulk Tray Dryer, and Triad

**7756100—Tray with Slide-Out Bottom.** Stainless steel. **Not for use with Small Tray Dryer.** 12.0" w x 14.0" d. Shipping weight 5.0 lbs. (2.3 kg)



**7369600—Microwell Plate Holder.** Holds a 127.8 x 85.5 mm microplate. 5.7" w x 8.0" d x 0.7" h. One each. Aluminum. **Microplate not included.** Shipping weight 2.0 lbs. (0.9 kg)



**Rectangular Test Tube & Microcentrifuge Tube Racks.** 304 stainless steel. 12.0" w x 10.0" d x 2.5" h. Holds up to 70 each 12-13 mm test tubes or microcentrifuge tubes. **Tubes not included.**

**7370000—Test Tube Rack.** Shipping weight 2.8 lbs. (1.3 kg)

**7370001—Microcentrifuge Tube Rack.** Includes adapter plate for microcentrifuge tubes. Shipping weight 3.5 lbs. (1.6 kg)



**7439300—Microcentrifuge Tube Holder.** 3.0" w x 2.25" d x 1.0" h. Holds twelve 1.7 ml tubes. Anodized aluminum. One each. May be placed on any flat freeze drying surface. **Tubes not included.** Shipping weight 0.6 lb. (0.3 kg)



**7516200—Support Grid.** Stainless steel tray with removable plastic grid holds 144 containers of 12 mm dia. 7.0" w x 7.0" d. Shipping weight 2.5 lbs. (1.1 kg)

# Accessories

## Manifolds

**Manifolds.** 304 stainless steel with neoprene valves and molded plastic knobs. For 1/2" and 3/4" flask adapters.



**7201700—4-Port.** 8.6" w x 8.6" d x 11.2" h.  
Shipping weight 5.5 lbs. (2.5 kg)



**7201800—8-Port.** 8.6" w x 8.6" d x 15.7" h.  
Shipping weight 8.5 lbs. (3.9 kg)



**7522300—20-Port.** 27.4" w x 8.7" d x 11.0" h.  
Shipping weight 15 lbs. (7 kg)



**7201900—14-Port.** For Console models only.  
24.3" w x 10.6" d x 10.3" h.  
Shipping weight 21 lbs. (10 kg)



**7202000—28-Port.** For Console models only.  
24.3" w x 10.6" d x 18.3" h.  
Shipping weight 38 lbs. (17 kg)

# End-Zone™ Starter Kit & Add-On Kit

## Specifications & Ordering Information

**End-Zone** provides an alternative to visual inspection, eliminating the guesswork of knowing when lyophilization is complete. End-Zone alerts when sample moisture levels are below your specified value. Lyo-Works shows the end point alert on your freeze dryer's display and, if programmed, sends an email.

End-Zone Starter Kit 780901000 connected to FreeZone 4.5 Liter Benchtop Freeze Dryer 700401000 with 8-Port Manifold 7201800



**Starter Kit** has all the components needed to detect end point in one 600 ml or larger Fast-Freeze Flask:

- 26-pin cable for connecting End-Zone sensors to any FreeZone Freeze Dryer with Lyo-Works (freeze dryer sold separately)
- Holder mountable to the freeze dryer with space for up to six sensors
- Two vacuum sensors—a sample sensor for placement in the Fast-Freeze Flask and a reference sensor for placement in an empty valve port
- Valve adapter to facilitate connection of sensor to empty valve port
- Fast-Freeze Flask top with opening for insertion of the vacuum sensor (flask bottom is required)
- Up to four End-Zone Add-On Kits may be used at once

**780901000—End-Zone Starter Kit.** Includes powder-coated steel 6-position sensor holder, two vacuum sensors with communication cables, one valve adapter and one top for 600 ml or larger Fast-Freeze Flask. FreeZone Freeze Dryer with Lyo-Works OS, drying chamber or manifold, vacuum pump, 600 ml or larger Fast-Freeze Flask and Adapter are required for operation. **Not for use with Stoppering or Bulk Tray Dryers.** Shipping weight 3.0 lbs. (1.4 kg)

**780902000—End-Zone Add-On Kit.** Includes one vacuum sensor with communication cable, one valve adapter, one top for 600 ml or larger Fast-Freeze Flask. Requires 780901000. Shipping weight 3.0 lbs. (1.4 kg)

# Accessories

## Stoppering Chambers & Optional Lids



**7569000—Mini Stoppering Chamber.** For freeze drying samples in small bottles with split stoppers. Attaches to valve on any Labconco Drying Chamber or Manifold, to stopper bottles under vacuum. For up to forty-eight 2 ml bottles or eight 30 ml serum bottles. Stoppering travel distance is 0.5". Using zero, one or two of the spacers provided, the chamber can accommodate bottles that range in height from 1.38" to 2.48". Includes one tray, two spacers; 1" OD and 1/2" ID PVC tubing, 36" long; two hose clamps; and 1/2" dia. stainless steel connector. Shipping weight 15 lbs. (7 kg)

@ **WARNING:** Cancer—P65Warnings.ca.gov (California only)



**Clear Stoppering Chambers.** For an economical means of stoppering serum bottles under vacuum. Includes acrylic chamber and 3/4" thick attachment port base, stainless steel top plate, two 10" dia. aluminum shelves and low voltage transformer with variable heat control to +40° C (+104° F). 3/4" thick base has neoprene gasket and pass-through for electrical transformer cord. When lyophilization is complete, handle is turned to collapse top plate onto shelves, stoppering samples under vacuum. For use with small containers from 2 - 5 ml.

**7868020—115V, 50/60 Hz.** Shipping weight 30 lbs. (14 kg)

**7868030—230V, 50/60 Hz.** Shipping weight 30 lbs. (14 kg)



7207500

**Optional Lids.** For added protection against solvents and corrosives that attack acrylic, such as acetonitrile. Lid has 3" dia. hole for direct connection of a chamber or manifold to the collector. Includes cap for closing port when not in use.

**7207400—Stainless steel.** 14.0" dia.

Not for use with 2.5 L models.

Shipping weight 5.0 lbs. (2.3 kg)

**7207500—Stainless steel.** 10.1" dia.

For use with 2.5 L models only.

Shipping weight 5.0 lbs. (2.3 kg)

**7208200—Acrylic.** 14.0" dia.

Shipping weight 5.0 lbs. (2.3 kg)

**7439902—Glass.** 14.0" dia.

For use with 16-Port Drying Chamber and 6, 12 and 18 L models only.

Shipping weight 3.5 lbs. (1.6 kg)

**7439900—Glass.** 9.25" dia.

For use with 12-Port Drying Chamber.

Shipping weight 2.0 lbs. (0.9 kg)



7208200

# Accessories

## Secondary Traps



**Secondary Traps.** Fits between the freeze dryer and vacuum pump to prevent particulate contaminants from entering the pump.

@ **WARNING:** Cancer—P65Warnings.ca.gov (California only)

**7386500—HEPA Filter Trap.**

Shipping weight 15 lbs. (7 kg)

**Chemical Trap.** Connect between the freeze dryer and vacuum pump for additional pump protection.



**7815302—Clear Polypropylene Canister with stand.**

For use with any trap insert listed below. Shipping weight 5.0 lbs. (2 kg)

**Trap Inserts.** Shipping weight 4.0 lb. (1.8 kg)

**7814800—Acid.** Contains acid neutralizing desiccant which changes from off-white to purple when exhausted.

**7814900—Moisture.** Contains desiccant which changes from blue to pink when exhausted.

**7815000—Radioisotope.** Contains activated carbon.

**7815200—Solvent.** Contains molecular sieve.

**7995600—Ammonia.** Contains impregnated carbon.



**7394800—Liquid Nitrogen Secondary Trap.** For samples with ultra low eutectic points. Liquid nitrogen cools to  $-196^{\circ}\text{C}$  ( $-321^{\circ}\text{F}$ ) to trap contaminants. Powder-coated steel.

3/4" OD inlet and outlet connections. Includes insulated filler tube for connection of user-supplied liquid nitrogen tank, 20" wire-reinforced PVC tubing, 3/4" ID and two clamps.

**Liquid nitrogen is required (not included).** Shipping weight 24 lbs. (11 kg)

@ **WARNING:** Cancer—P65Warnings.ca.gov (California only)



**7538000—Dry Ice Secondary Trap.** For samples with low eutectic points. Dry ice and solvent cools the insulated well to approximately  $-75^{\circ}\text{C}$  ( $-103^{\circ}\text{F}$ ). When used alone, it may serve as an inexpensive collector. Connect a sample directly to the incoming port and add a vacuum pump to the outgoing port. Ice trapping capacity of 900 ml. Includes 3.1 L well and 3/4" OD inlet and outlet connections. **Dry ice and solvent, such as alcohol, are required (not included).** Shipping weight 9.0 lbs. (4.1 kg)

# Accessories

## Fast-Freeze™ Flasks

**Fast-Freeze Flasks.** Flasks should be filled to no more than one-third of their volume. Includes silicone top, glass bottom and a supply of filter paper. Amber-colored glass helps to protect light-sensitive samples and is available in four sizes. Clear available in 10 sizes.

**Adapters are required (not included).**

- Compatible with all major brands of laboratory freeze dry apparatus
- Variety of adapters available
- No complex filter retainer is necessary. One hundred filters supplied with each flask
- High strength borosilicate glass bottom to withstand extreme temperatures and high vacuum
- Entire flask may be autoclaved



| Flask Size | Catalog Number |                | Flask Top Diameter* | Flask Bottom Dimensions (h x d) | Shipping Weight   |
|------------|----------------|----------------|---------------------|---------------------------------|-------------------|
|            | Clear Flask    | Amber Flask    |                     |                                 |                   |
| 40 ml      | <b>7540000</b> |                | 1/2"                | 76 x 34 mm                      | 1.0 lb. (0.5 kg)  |
| 80 ml      | <b>7540200</b> |                | 1/2"                | 115 x 34 mm                     | 1.0 lb. (0.5 kg)  |
| 120 ml     | <b>7540300</b> |                | 3/4"                | 68 x 59.2 mm                    | 1.0 lb. (0.5 kg)  |
| 150 ml     | <b>7540400</b> |                | 3/4"                | 85 x 59.2 mm                    | 1.0 lb. (0.5 kg)  |
| 300 ml     | <b>7540600</b> | <b>7540601</b> | 3/4"                | 145 x 59.2 mm                   | 1.5 lbs. (0.7 kg) |
| 600 ml     | <b>7540800</b> | <b>7540801</b> | 3/4"                | 135 x 90.2 mm                   | 2.0 lbs. (0.9 kg) |
| 750 ml     | <b>7541100</b> | <b>7541101</b> | 3/4"                | 162 x 90.2 mm                   | 2.0 lbs. (0.9 kg) |
| 900 ml     | <b>7540900</b> | <b>7540901</b> | 3/4"                | 190 x 90.2 mm                   | 2.2 lbs. (1.0 kg) |
| 1200 ml    | <b>7541000</b> |                | 3/4"                | 240 x 90.2 mm                   | 2.2 lbs. (1.0 kg) |
| 2000 ml    | <b>7541200</b> |                | 3/4"                | 380 x 90.2 mm                   | 2.2 lbs. (1.0 kg) |



**7544810—Replacement Filter Paper.** 1000 (10 packages of 100). 1" dia. (28mm). Pore size 17 microns. Shipping weight 0.5 lb. (0.2 kg)

\*Valves on the freeze dryers, drying chambers and manifolds fit adapters with both 1/2" and 3/4" dia.

# Accessories

## Fast-Freeze™ Adapters

**Flask Adapters.** Connect Fast-Freeze Flasks to the valve ports on your drying chamber or manifold. Borosilicate glass or stainless steel adapters in 1/2" and 3/4" dia. **Required for operation.** Shipping weight 0.5 lb. (0.2 kg)



| Straight Stainless Steel |              |
|--------------------------|--------------|
| <b>7547000</b>           | 1/2" to 1/2" |
| <b>7547200</b>           | 3/4" to 3/4" |

| Straight Borosilicate Glass |                               |
|-----------------------------|-------------------------------|
| <b>7545000</b>              | 1/2" to 1/2"                  |
| <b>7545200</b>              | 3/4" to 3/4"                  |
| <b>7545400</b>              | 1/2" to 3/4" and 3/4" to 1/2" |

| 45° Bend Stainless Steel |              |
|--------------------------|--------------|
| <b>7547400</b>           | 1/2" to 1/2" |
| <b>7547600</b>           | 3/4" to 3/4" |

| 45° Bend Borosilicate Glass |                               |
|-----------------------------|-------------------------------|
| <b>7545600</b>              | 1/2" to 1/2"                  |
| <b>7545800</b>              | 3/4" to 3/4"                  |
| <b>7546000</b>              | 1/2" to 3/4" and 3/4" to 1/2" |



**7543800—Stainless Steel Adapter 29/32.** Allows connection of 29/32 round bottom flasks to valve ports on drying chambers and manifolds. Shipping weight 1.0 lb. (0.5 kg)



**7543900—Aseptic Adapter.** For use with Fast-Freeze Flasks in place of an adapter. Allows a flask to be filled under sterile conditions prior to lyophilization and transported in a closed position before and after. 3/4" OD for connection between a Fast-Freeze Flask and a valve port on any Labconco Drying Chamber or Manifold. Corrosion-resistant PTFE. 1.75" w x 0.75" d x 3.2" h. Shipping weight 0.5 lb. (0.2 kg)



**7540100—FreeZone Freeze Dry Glassware Kit.** Includes one each flask and adapters: 7540800 (600 ml) and 7540600 (300 ml), 7547200 (3/4" straight stainless steel) and 7547600 (3/4" 45° bend stainless steel). Shipping weight 5.0 lbs. (2.0 kg)

# Accessories

## Fast-Freeze™



**Tube Holders.** Stainless steel. Hold tubes in an upright position inside a Fast-Freeze Flask or Lymph-Lock Flask. **Tubes and glassware are not included.** Shipping weight 1.0 lb. (0.5 kg)

**7442600—Holds 15 tubes.** For 12 or 13 mm dia tubes. 144 h x 83.3 d mm. For flasks 900 ml or larger.

**7442700—Holds 30 tubes.** For 0.5, 1.0 and 2.0 ml microcentrifuge tubes. 139 h x 83.3 d mm. For flasks 900 ml or larger.

**7379400—Holds 9 tubes.** For 15 ml tubes. 170 h x 88.9 d mm. For flasks 1200 ml or larger.

**7379500—Holds 4 tubes.** For 50 ml tubes. 30 h x 116 d mm. For flasks 1200 ml or larger.



**7369700—Slant Freeze Flask Holder.** Slant freezing increases the surface area to volume ratio, maximizing freeze drying rate. Angle is manually adjustable from 5° to 35°. Powder-coated steel. 10.4" w x 4.5" d x 6.7" h. Shipping weight 3.5 lbs. (1.6 kg)



**Flask Holders.** Help to prevent sample meltback. Hold Fast-Freeze Flask containing the sample inside a larger Fast-Freeze Flask so that sample is insulated. Stainless steel. **Flasks not included.** Shipping weight 1.0 lb. (0.5 kg)

**7543300—Small.** For use with inner flasks 7540000 (40 ml), 7540200 (80 ml) and outer flasks 7541100 (750 ml), 7541101 (750 ml), 7540900 (900 ml), 7540901 (900 ml).

**7543400—Large.** For use with inner flasks 7540300 (120 ml), 7542400 (150 ml), 7540600 (300 ml), 7540601 (300 ml) and outer flask 7541000 (1200 ml).

# Accessories

## Serum Bottles & Ampules

**Serum Bottles.** Specifically designed for lyophilization with uniform thin walls to even freezing and drying. Serum Bottles, Stoppers and Seals supplied in packages of 100.




| Size   | 20 mm Corkage | 13 mm Corkage | Split Stoppers | Aluminum Seals |
|--------|---------------|---------------|----------------|----------------|
| 2 ml   |               | 7575010       | 7576010        | 7577010        |
| 3 ml   |               | 7575210       | 7576010        | 7577010        |
| 5 ml   | 7573010       |               | 7576210        | 7577110        |
| 10 ml  | 7573210       |               | 7576210        | 7577110        |
| 20 ml  | 7573410       |               | 7576210        | 7577110        |
| 30 ml  | 7573610       |               | 7576210        | 7577110        |
| 50 ml  | 7573810       |               | 7576210        | 7577110        |
| 100 ml | 7574010       |               | 7576210        | 7577110        |
| 125 ml | 7574210       |               | 7576210        | 7577110        |



**Seal Crimpers.** Secure tear-away Aluminum Seals. Shipping weight 3.0 lbs. (1.4 kg)

7578000—For 13 mm corkage.

7578100—For 20 mm corkage.

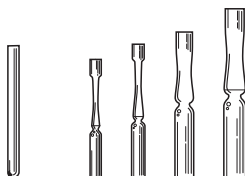

**Ampules.**  
Borosilicate glass.  
Package of 100.

| Size  | Straight | Pre-Scored Flat Bottom |
|-------|----------|------------------------|
| 1 ml  | 7572710  | 7570610                |
| 2 ml  |          | 7570810                |
| 5 ml  |          | 7571010                |
| 10 ml |          | 7571210                |



**Ampule Pods.** Attach to a single drying chamber or manifold valve port. Polypropylene base with 15 silicone tubing ports, 15 plugs and two clamps. Include one straight stainless steel adapter for connection to a drying accessory valve port. 3.8" w x 3.8" d x 7.5" h. Shipping weight 2.5 lbs. (1.1 kg)

7205100—1/4" ID tubing port dia. For use with 7572710 and 7571010.

7205101—3/8" ID tubing port dia. For use with 7570610, 7570810 and 7571210.



**7593401—Ampule Valve Adapters.** Connect ampules to valve ports.  
Package of 10. Shipping weight 0.5 lb. (0.2 kg)



**7762700—Three Way Adapter.** Permits attachment of three ampules to a single valve port. Ampules are attached using Ampule Valve Adapters 7593401 or 1/8" surgical tubing (**valve adapters and tubing not provided**). Cavity in adapter body can be filled with cotton fiber media to help prevent contamination between samples. Shipping weight 0.2 lb. (0.1 kg)



**7578500—Oxygen/Natural Gas Sealing Torch.** Specifically designed for flame sealing freeze dry ampules. Seals all types of heat-resistant glass. Connects to natural gas, butane or propane and oxygen with 1/4" ID hose connectors. Shipping weight 3.0 lbs. (1.4 kg)

# FreeZone® Benchtop Shell Freezers

## Specifications & Ordering Information



### All models feature:

- Stainless steel and powder-coated steel exterior
- 1/2 hp HCFC/CFC-free refrigeration system to cool the heat transfer solution to -40° C (-40° F)
- Stainless steel bath, 5.5" w x 12.0" d x 6.5" h, surrounded by urethane foam insulation with white high-density polyethylene cover. Two chain-driven rollers, powered by a sparkless induction motor, rotate flasks up to 1200 ml. Approximately two liters of transfer solution are required (not included)
- LED temperature "wave" for display of bath temperature from -20° to -42° C (-4° to -43.6° F). LEDs illuminate amber as the system cools. When system reaches operating temperature, bottom LED illuminates green
- Retractable drain hose
- Fast-Freeze Flask capacity:
  - Two: 40 ml, 80 ml, 120 ml and 150 ml
  - One: 300 ml, 600 ml, 750 ml, 900 ml and 1200 ml
- Overall dimensions: 15.1" w x 22.1" d x 19.7" h
- Shipping weight 96 lbs. (44 kg)

### FreeZone Benchtop Shell Freezers.

For pre-freezing samples. Stainless steel bath has two rollers that rotate Fast-Freeze Flasks in a low temperature heat transfer solution, such as methanol or ethanol, so that a uniform thin coating of sample freezes to the interior surface of the container. Shell freezing increases the surface area to volume ratio, which facilitates efficient vapor flow from sample to collector during lyophilization.

### Models conform to the following standards:

- UL Standard 61010-1 (115/230V, 60 Hz models)
- CAN/CSA C22.2 No. 61010-1 (115/230V, 60 Hz models)
- CE Conformity marking (230V, 50 Hz models)

### Option includes:

- Region-specific electrical configurations

### Standard Refrigerant

| Catalog Number | Electrical Requirements | Plug Type                |
|----------------|-------------------------|--------------------------|
| 7949020        | 115V, 60 Hz, 6A         | North America, 115V, 15A |
| 7949040        | 230V, 60 Hz, 3A         | North America, 230V      |
| 7949030        | 230V, 50 Hz, 3A         | Schuko                   |
| 7949060        | 230V, 50 Hz, 3A         | British (UK)             |
| 7949070        | 230V, 50 Hz, 3A         | China/Australia          |

### Natural Refrigerant

| Catalog Number | Electrical Requirements | Plug Type                                                                   |
|----------------|-------------------------|-----------------------------------------------------------------------------|
| 794905090      | 230V, 50 Hz, 3A         | Multi-Plug: Schuko, British (UK), China/Australia, and India plugs included |



# Accessories

## Carts, Shelves & Gauges



**Digital Electronic Vacuum Gauges.** LCD displays vacuum measurement in mBar, Torr or Pa. Range is 0.02 to 5 mBar. Stainless steel and powder-coated steel casing. 8.1" w x 4.1" d x 6.4" h. Include a cord with 3/8" to 1/2" OD adapter for connection to either vacuum tubing or a valve on a drying accessory. May be used as an accessory by the Dry Ice Benchtop Freeze Dryer or other freeze dryers. For use with Dry Ice Benchtop Freeze Dryer.

Shipping weight 10 lbs. (5 kg)

**7395000—115V, 60 Hz, 0.4A.**

**7395001\*—230V, 50 Hz, 0.2A.**

\*International electrical configuration



**7368700—FreeZone Cart.** Powder-coated steel. 5.0" dia. casters. Bottom shelf accommodates a vacuum pump. 23.3" between shelves. Load capacity of 450 lbs. (204 kg). 28.0" w x 28.0" d x 31.3" h. Shipping weight 120 lbs. (54 kg)



**7426400—Support Shelves.** Attach easily to Benchtop Freeze Dryers to support large flasks attached to the drying chamber or manifold. Shelf height is adjustable. One pair. Powder-coated steel. Each shelf is 7.4" w x 12.0" d. Shipping weight 2.0 lbs. (0.9 kg)

# FreeZone® Freeze Dryers

## Specifications

### All models feature:

- ❖ Full color, 5" capacitive touch screen with Lyo-Works OS:
  - Real time display of collector temperature, vacuum level and optional End-Zone end point detection
  - Automatic and manual start-up modes for collector cool down and vacuum pull down
  - Vacuum set point, alerts, diagnostic tests, user manual, defrost control and optional drying accessory control
  - Data logging stores and displays data in table or graph format, and may be exported via USB
  - Settings include language (English, Spanish, French, German, Italian, Portuguese, Chinese, and Japanese), time, temperature (°C or °F) and vacuum (mBar, Torr or Pa)
  - Alerts can be emailed to user-entered addresses via Ethernet connection
  - Diagnostics and QR Resources
  - Maintenance alert settings for vacuum pump and refrigeration
  - Security lock-out
- ❖ USB port and Ethernet connection
- End-Zone compatible (End-Zone Starter Kit sold separately. (See page 26))
- Vacuum control valve and break valve
- ❖ Patented\* vacuum line moisture sensor
- Stainless steel and powder-coated steel exterior
- Quick-disconnect drain hose fitting with hose connector. 18" of 1/4" ID silicone tubing and clamp included
- Collectors reaching -50°C include acrylic lid. Collectors reaching -84°C or -105°C include stainless steel lid
- Hot gas defrost with auto shut off at +60°C (+140°F)

### Benchtop models also feature:

- Upright stainless steel collector coil capable of holding 2.5, 4.5 or 8 liters of ice before defrosting. Collector lid, 3/4" thick, has neoprene gasket and 3" dia. port to connect drying accessories (sold separately)
- Electrical receptacle for vacuum pump connection
- 3/4" OD vacuum connection, 50" of 3/4" ID vacuum hose and two clamps

### Console models also feature:

- Upright stainless steel collector coil capable of holding 6, 12 or 18 liters of ice before defrosting. Collector lid, 3/4" thick, has neoprene gasket
- Console cabinet with interior space for vacuum pump (sold separately) and two-piece, removable front panel with viewing window
- Clear acrylic chamber lid, 3/4" thick, with neoprene gasket
- Exterior electrical receptacle for vacuum pump connection

- 3/4" OD vacuum connection, 36" of 3/4" ID vacuum hose and two clamps

### Standards conformance & regulations:

- UL Standard 61010-1 (115/230V, 60 Hz models)
  - CAN/CSA C22.2 No. 61010-1 (115/230V, 60 Hz models)
  - CE Conformity marking (230V, 50 Hz models)
- @ **WARNING:** Cancer — P65Warnings.ca.gov (California only)

### Benchtop options include:

- **For 2.5 and 4.5 L:** Single 1/3 hp HCFC/CFC-free refrigeration system to cool collector to ~-50°C (~-58°F), two 1/3 hp HCFC/CFC-free refrigeration systems to cool collector to ~-84°C (~-119°F) or **for 4.5 L only:** two 1/2 hp HCFC/CFC-free refrigeration systems to cool collector to ~-105°C (~-157°F) at 60 Hz or ~-102°C (~-152°F) at 50 Hz. **For 8 L:** Single 3/4 hp HCFC/CFC-free refrigeration system to cool collector to ~-50°C (~-58°F)
- PTFE-coated collector coil and chamber for corrosive compounds
- Region-specific electrical plug configurations

### Console options include:

- Single 3/4 hp (6 L), 1 hp (12 L) or 1-1/2 hp (18 L) HCFC/CFC-free refrigeration system to cool collector to ~-50°C (~-58°F) or two 1/3 hp (6 L) or 3/4 hp (12 L) HCFC/CFC-free refrigeration systems to cool collector to ~-84°C (~-119°F)
- Purge valve - temporarily isolates vacuum from the system to optimize performance before or after freeze drying
- PTFE-coated collector coil and chamber for corrosive compounds
- Built-in shell freezer with stainless steel bath, 5.5" w x 12.5" d x 7.5" h, that includes 1/2 hp HCFC/CFC-free refrigeration system to cool heat transfer solution to -40°C (-40°F), 3/4" thick high-density polyethylene lid, neoprene gasket and quick-disconnect drain hose fitting and 18" of 1/4" ID silicone tubing with clamp (6 L -50°C only)
- Built-in stainless steel mini vacuum drying chamber, 5.1" w x 13.0" d x 2.5" h, that includes 50-watt heater with microprocessor control to +60°C (+140°F), 3/4" thick clear acrylic lid and neoprene gasket. (Chamber is not cooled.)
- Region-specific electrical plug configurations

### All models require (not included):

- Vacuum pump with a displacement of at least 98 L/min, 0.02 mBar ultimate pressure and fitting for 3/4" ID vacuum hose. See pages 19-20
- Drying accessory. See pages 14-15, 17-18, 22, 24-25
- Freeze dry glassware if not bulk freeze drying. See pages 29-32

❖ Labconco exclusive feature  
\*U.S. Patent No. 6,968,631



# FreeZone® Triad® Freeze Dryers

## Specifications

### All models feature:

- \* Stainless steel collector coil capable of removing 1.84 L of water in 24 hours and holding 2.5 L of ice before defrosting\*
- \* Two 1/3 hp HCFC/CFC-free refrigeration compressors cool the collector to  $\approx -85^{\circ}\text{C}$  ( $\approx -121^{\circ}\text{F}$ ) and work in concert with a 1000-watt heater to cool and heat fluid medium circulating through channels in the shelf. Not for use with samples containing methanol or ethanol
- Processing shelf, 12.4" w x 14.5" d. Spacing between shelf and top of chamber accommodates containers with stoppers up to 148 mm maximum height
- One probe for monitoring temperature of sample in chamber
- Temperature ramp and hold programming from  $\approx -55^{\circ}\text{C}$  to  $\approx +50^{\circ}\text{C}$  ( $\approx -67^{\circ}\text{F}$  to  $\approx +122^{\circ}\text{F}$ ) and  $\approx -75^{\circ}\text{C}$  ( $\approx -103^{\circ}\text{F}$ ) during pre-freezing. Controls temperature to  $\pm 1^{\circ}\text{C}$  of set point. Memory to store 30 programs with 16 segments each including number of steps, ramp rate, shelf temperature, hold time and vacuum level
- \* Full color, 5" capacitive touch screen with Lyo-Works OS:
  - Real time display of collector temperature and vacuum level
  - Creation and storage for up to 30 programs with up to 16 steps each. Ability to set ramp rate, shelf temperature, shelf temperature time and vacuum set point per step
  - Manual start-up buttons for collector and shelf temperature control and vacuum pull down
  - Alerts, diagnostic tests, user manual and defrost control
  - Data logging stores and displays data in table or graph format, and may be exported via USB or Ethernet
  - Settings include language (English, Spanish, French, German, Italian, Portuguese, Chinese, and Japanese), time, temperature ( $^{\circ}\text{C}$  or  $^{\circ}\text{F}$ ) and vacuum (mBar, Torr or Pa)
  - Alerts can be emailed to user-entered addresses via Ethernet connection
  - Vacuum pump maintenance alert
  - Security lock-out
  - Diagnostics
  - QR Resources
- \* Pneumatic stoppering mechanism with diaphragm and plate.  
**No compressed gas is required**
- Stoppering control knob, vacuum release control knob and 1/8" OD user-supplied sterile/inert gas backfill port
- Four left side-mounted neoprene valves with molded plastic knobs and fittings that accommodate 1/2" and 3/4" flask adapters
- \* USB port and Ethernet connection
- Electrical receptacle for vacuum pump
- Hot gas defrost with stainless steel drain pan, silicone drain hose and auto shut off at  $+60^{\circ}\text{C}$  ( $+140^{\circ}\text{F}$ )
- Vacuum control valve
- \* Vacuum break valve
- Stainless steel and powder-coated steel exterior
- Clear acrylic door, 1" thick, with silicone gasket and latch
- 3/4" OD vacuum connection. 36" long, 3/4" ID vacuum hose with two clamps

### Standards conformance & regulations:

- UL Standard 61010-1 (115/230V, 60 Hz models)
- CAN/CSA C22.2 No. 61010-1 (115/230V, 60 Hz models)
- CE Conformity marking (230V, 50 Hz models)

@ **WARNING:** Cancer — P65Warnings.ca.gov (California only)

### Options include:

- Region-specific electrical plug configurations

### All models require (not included):

- Vacuum pump with a displacement of at least 98 L/min, 0.02 mBar ultimate pressure and fitting for 3/4" ID vacuum hose. See pages 19-20
- Freeze dry glassware if not bulk freeze drying. See pages 29-32

\* Freeze drying rate will be lower for samples other than plain water.  
For optimum performance, room temperature should be  $21^{\circ}\text{C}$  ( $70^{\circ}\text{F}$ ) or colder.

\* Labconco exclusive feature



# FreeZone® Stoppering Tray Dryers

## Specifications

### All models feature:

- Three adjustable shelves, each 12.0" w x 16.5" d with 196 square inches of area. Spacing between the shelves accommodates containers with stoppers up to a maximum height of 75 mm and minimum height of 38 mm when three shelves are used, a maximum height of 107 mm and minimum height of 56 mm when two shelves are used, and a maximum height of 196 mm and minimum height of 114 mm when one shelf is used
- 1 hp HCFC/CFC-free refrigeration system and 1000-watt electric heater for cooling and heating fluid medium circulating through channels in the shelves from -40° C to +40° C (-40° F to +104° F). Controls circulating fluid temperature to  $\pm 1^\circ$  C of set point
- Three probes for monitoring sample or shelf temperature
- Temperature ramp and hold programming from -40° C to +40° C (-40° F to +104° F) and memory to store 30 programs with 16 segments each, including number of steps, ramp rate, shelf temperature, hold time and vacuum level
- \* System programming that is displayed on the FreeZone 6, 12 or 18 Liter Freeze Dryer (freeze dryer sold separately or as one catalog number): Manual or Program mode; probe temperatures in °C or °F; vacuum in mBar, Torr or Pa; segment of the program that is in progress; time remaining in present segment; alarms for power failure; supply line voltage outside of limits; system temperature variation  $> \pm 2^\circ$  C, vacuum variation progress bar, diagnostics, and QR resources  $> 0.500$  mBar; failure to reach system set point temperature; and program completed
- \* Data logging that displays and stores data in table or graph on the FreeZone 6, 12 or 18 Liter Freeze Dryer (freeze dryer sold separately or as one catalog number). Export data via USB or Ethernet
- Stoppering mechanism that inflates a diaphragm via atmospheric pressure. **No compressed gas is required**
- Stoppering Raise/Lower control knob, vacuum release Open/Close control knob and 1/8" OD backfill port for user-supplied air or gas
- Stainless steel and powder-coated steel exterior
- Acrylic door, 1" thick, with silicone gasket and latch
- 3" dia. outlet, located on the bottom of the tray dryer cabinet, for connection to a FreeZone 6, 12 or 18 Liter Console Freeze Dryer (freeze dryer is required and not included)
- Integral support stand
- Overall dimensions (tray dryer only): 32.4" w x 24.6" d x 27.1" h

### Standards conformance & regulations:

- UL Standard 61010-1 (115/230V, 60 Hz models)
- CAN/CSA C22.2 No. 61010-1 (115/230V, 60 Hz models)
- CE Conformity marking (230V, 50 Hz models)

### Options include:

- 6-Port Tray Dryer Manifold, 30.0" w x 16.5" d x 7.0" h. Available only when ordering the FreeZone Stoppering Tray Dryer as a separate component. See page 21
- Isolation Valve. Available only when ordering the FreeZone Stoppering Tray Dryer as a separate component. See page 21
- Region-specific electrical plug configurations

### All models require (not included):

- FreeZone 6, 12 or 18 Liter Console Freeze Dryer with Lyo-Works OS. FreeZone Stoppering Tray Dryer and Freeze Dryer may be ordered as separate components or together as one catalog number on selected models of FreeZone 6, 12 or 18 Liter Console Freeze Dryers. See page 11
- Vacuum pump with a displacement of at least 163 L/min, 0.02 mBar ultimate pressure and fitting for 3/4" ID vacuum hose. See pages 19-20
- Freeze dry glassware if not bulk freeze drying. See pages 29-32

\* Labconco exclusive feature



# FreeZone® Bulk & Small Tray Dryers

## Specifications

### All models feature:

- ❖ 200-watt heater on each shelf for heating to +60° C (+140° F). Maintains system temperature within 3° C of set point. (Shelves are not cooled by the Tray Dryer. Samples must be pre-frozen)
- Three sensor probes for monitoring sample temperature. Connections for two additional sensors are provided (additional sensors sold separately)
- ❖ Data logging, when enabled, stores and displays data in table or graph on the FreeZone Freeze Dryer. (Freeze dryer sold separately. Export data via USB or Ethernet)
- Vacuum release valve with knob
- Stainless steel and powder-coated steel exterior
- Acrylic door, 1" thick, with silicone gasket and latch

### Bulk Tray Dryers also feature:

- Three shelves, each 12.7" w x 16.6" d to provide 210 square inches of area. Chamber may accommodate two additional shelves (additional shelves sold separately)
- ❖ System programming displayed on the FreeZone 6, 12 or 18 Liter Freeze Dryer (freeze dryer sold separately): Manual or Program mode; probe temperatures in °C or °F; vacuum in mBar, Torr or Pa; segment of the program in progress; time remaining in present segment; alarms; shelf temperature variation  $>\pm 3^{\circ}\text{C}$ ; failure to reach shelf set point temperature  $\pm 3^{\circ}\text{C}$ ; and program completed
- 3" dia. outlet, located on the bottom of the tray dryer cabinet, for connection to a FreeZone 6, 12 or 18 Liter Console Freeze Dryer (freeze dryer sold separately)
- ❖ Integral support stand with attachment port for connection to FreeZone 6, 12 or 18 Liter Console Freeze Dryer (freeze dryer sold separately)
- Overall dimensions: 32.2" w x 21.5" d x 27.1" h

### Small Tray Dryers also feature:

- Three removable shelves with shelf sensors, each 10.6" w x 12.0" d to provide 127.2 square inches of area. Chamber can accommodate two additional shelves (additional shelves sold separately)
- ❖ System programming is displayed on any FreeZone 2.5, 4.5, 6, 8, 12 or 18 Liter Freeze Dryer (Freeze dryer sold separately. Manual or Program mode; probe temperatures in °C or °F; vacuum in mBar, Torr or Pa; segment of the program in progress; time remaining in present segment; end of program; alarms; supply line voltage outside of limits; system temperature variation  $>\pm 2^{\circ}\text{C}$ , vacuum variation  $>0.500\text{ mBar}$ ; failure to reach system set point temperature; and program completed
- Overall dimensions: 15.3" w x 18.5" d x 18.1" h

### Standards conformance & regulations:

- ❖ Labconco exclusive feature

- UL Standard 61010-1 (115/230V, 60 Hz models)
- CAN/CSA C22.2 No. 61010-1 (115/230V, 60 Hz models)
- CE Conformity marking (230V, 50 Hz models)

### Bulk Tray Dryer options include:

- 6-Port Tray Dryer Manifold, 30.0" w x 16.5" d x 7.0" h
- Isolation Valve

### Bulk Tray Dryer models require (not included):

- FreeZone 6, 12 or 18 Liter Console Freeze Dryer with Lyo-Works OS. See page 10-11
- Vacuum pump with a displacement of at least 163 L/min, 0.007 mBar ultimate pressure and fitting for 3/4" ID vacuum hose. See pages 19-20
- Freeze dry glassware if not bulk freeze drying. See pages 29-320

### Small Tray Dryer models require (not included):

- FreeZone 2.5 L, 4.5 L, 6 L, 8 L, 12 L, 18 L Freeze Dryer with Lyo-Works OS
- Vacuum pump with a displacement of at least 98 L/min, 0.02 mBar ultimate pressure and fitting for 3/4" ID vacuum hose. See pages 19-20
- Freeze dry glassware if not bulk freeze drying. See pages 29-32



**Visit [labconco.com](https://labconco.com) to request information  
for all Labconco products.**

**LABCONCO CORPORATION**

8811 Prospect Avenue  
Kansas City, MO 64132

(800) 821-5525 | (816) 333-8811

[labconco.com](https://labconco.com)



©2026 LABCONCO CORPORATION  
Product design subject to change without notice.

3-40-9/22/26

